

Now- what's the best choice in imaging systems? : Molecular Dynamics.

Contributors

Molecular Dynamics (Firm)

Publication/Creation

Kemsing : Molecular Dynamics, 1994.

Persistent URL

<https://wellcomecollection.org/works/ehakksjs>

License and attribution

Conditions of use: it is possible this item is protected by copyright and/or related rights. You are free to use this item in any way that is permitted by the copyright and related rights legislation that applies to your use. For other uses you need to obtain permission from the rights-holder(s).



Wellcome Collection
183 Euston Road
London NW1 2BE UK
T +44 (0)20 7611 8722
E library@wellcomecollection.org
<https://wellcomecollection.org>

What's the best choice?

**Molecular
Dynamics**
Molecular Dynamics, Ltd.
4 Chaucer Business Park
Kerning, Sevenoaks
Kent, TN11 0PL

**Molecular
Dynamics**

Molecular Dynamics, Ltd.
4 Chaucer Business Park
Kerning, Sevenoaks
Kent, TN11 0PL

**Now –
what's the best choice
in imaging systems?**

...Molecular Dynamics.

Two new high performance imaging systems from Molecular Dynamics, the Personal Densitometer™ SI and PhosphorImager™ SI, represent a breakthrough in image scanning and analysis. These powerful, laser-based scanners bring unprecedented flexibility and ease-of-use to phosphor screen imaging and densitometry. Both come with your choice of scanning and quantitative analysis software for either the Mac or PC-486.

The Mac and PC image analysis software includes:

- Area integration, including automatic peak finding
- Volume integration, including automatic band finding
- Multiple object drawing tools for precise measurements
- Convenient reporting of results
- Easy annotation of images

Molecular Dynamics' new Personal Densitometer SI offers the highest performance of any densitometer on the market. It features laser-based accuracy, high sensitivity and unmatched resolution.

- Linear range from 0 to 4 OD for high density readings
- Highest sensitivity for accurate reading of faint bands
- High-efficiency light collection for best quantitation
- Spatial resolution to 50 microns to resolve close bands
- Your choice of Mac or PC image analysis software





THE PERSONAL DENSITOMETER SI INCLUDES YOUR CHOICE OF MAC OR PC SOFTWARE.

Powerful Image Analysis Software for Filmless Autoradiography & Densitometry

NOW ALL THE
POWER OF THE
MACINTOSH IS
AVAILABLE ON
MOLECULAR
DYNAMICS'
IMAGING
SYSTEMS.

Mac

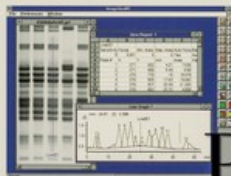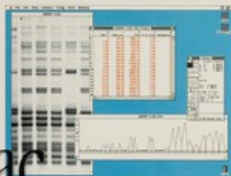

PC

WITH
THE NEW
IMAGEQUANT
SOFTWARE,
SCAN ONE
SAMPLE WHILE
ANALYZING
ANOTHER.



THE PHOSPHORIMAGER SI OFFERS THE FASTEST TECHNOLOGY FOR IMAGING RADIOACTIVE SAMPLES.

Molecular Dynamics

Innovative Instrument Systems for the Life Sciences

The PhosphorImager SI is a compact laser optical scanner that uses innovative storage phosphor screens to detect and quantitate radioactivity in gels, blots, TLC or tissues. With a linear dynamic range of four orders of magnitude, weak and strong signals require only a single exposure for complete visualization and analysis.

- Fast image generation without film
- Accurate quantitation without scintillation counting
- High throughput for individual labs or entire departments
- Your choice of Mac or PC image analysis software

**CALL TODAY FOR MORE DETAILS ON THE
BEST CHOICE IN IMAGE ANALYSIS.
0732-7625 65**

Fax: 0752-765422

Yes! I'm interested in Imaging Analysis.

For a complete information package on

the Molecular Dynamics Imaging Systems, please return the postage paid reply card.

☐ Personal Densitometer SI

☐ PhosphorImager SI

My platform is ☐ Mac ☐ PC (486)

☐ My request is urgent! Please expedite!

Name	_____
Address	_____
Post Office Box	_____
City	_____
State	_____
Zip	_____
Country	_____

Molecular Dynamics

Innovative Instrument Systems for the Life Sciences

© 1991 Molecular Dynamics. Molecular Dynamics and PhosphorImager are registered trademarks of Molecular Dynamics. All other trademarks and registered trademarks are the property of their respective owners.